

- Animal Welfare
- Flexibility
- Visibility
- Ergonomics





## Ferret Racks

### 1/Animal Welfare

Tecniplast racks for ferrets are compliant with US, UK and European Guidelines. Our two-levels system provides an area of 5400 cm<sup>2</sup> (5.8 ft<sup>2</sup>) with a smooth, warm floor and rounded corners. It includes important enrichment items, like raised shelves and an interlinking tube to connect two cages on the same rack. The tube can be open or closed to monitor the access of the animal.

### 2/Flexibility


It is easy to let several ferrets socialize, thanks to the removable transparent side panels and the possibility to connect two or more racks side by side. In just a few seconds your animals will have a bigger floor area for group housing and interactions.

### 3/Visibility

All the parts are designed to maximize the visibility of the animals inside the cage, with a large use of transparent polycarbonate, even at front door level.

### 4/Ergonomics

Both water bottle and feeder are external, for a very easy refilling. The door has a unique lockingsystem, for a safe and quick closure. The Noryl plastic cages are easy to wash, very light and can be frequently autoclaved at 121°C/250°F.

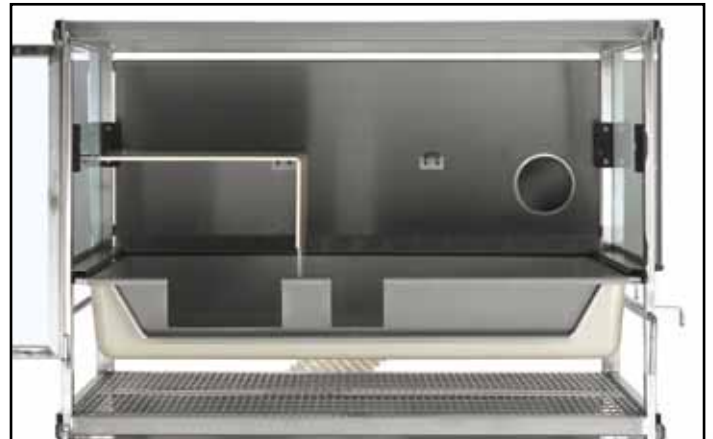

#### TECHNICAL SHEET

##### Cage Specification and Dimensions

###### 4541P008

"P" type solid base cage (compatible with tray Cat. N. 4501k389)  
Size (W x D x H): 870 x 740 x 130 mm / 34.25 x 29.13 x 5.118 inch  
Area (W x D x H): 797 x 679 x 122 mm / 31.38 x 26.73 x 4.803  
Weight: 3.8 Kg

##### Rack Specification and Dimensions

Complete Rack Unit Cat. N.: 4P02B700FT

"N. of Cages: 2  
Internal Height: 600 mm / 23.62 inch  
Overall Rack Dimensions (W x D x H):  
940 x 845 (+160) x 1700 mm  
37.01 x 33.27 (+6.299) x 66.93 inch



HEADQUARTERS **ITALY** | Tel. +39 0332 809711 • www.tecniplast.it • E-mail: tecnicom@tecniplast.it

SISTER COMPANIES

**AUSTRALIA/NEW ZEALAND** | Tel. + 61 2 8845 6500 • www.tecniplast.it • E-mail: info@tecniplast.com.au

**CHINA** | Tel. + 86 21 50810920 • www.tecniplast.cn • E-mail: tecniplastchina@tecniplast.it

**FRANCE** | Tél. 04 72 52 94 41 • www.tecniplast.fr • E-mail: info@tecniplast.fr

**GERMANY** | Tel. 08805 921320 • www.tecniplast.de • E-mail: info@tecniplast.de

**UNITED KINGDOM** | Tel. 0845 0504556 • www.tecniplast.it • E-mail: info@tecniplastuk.com

**JAPAN** | Tel. + 03 5770 5375 • www.tecniplast.it • E-mail: info@tecniplastjapan.co.jp

**USA/CANADA** | Toll Free: 877.669.2243 • www.tecniplastusa.com • E-mail: info@tecniplastusa.com

To find your local distributor, please visit [www.tecniplast.it](http://www.tecniplast.it)



ISO-9001 - Cert. n° 0875



ISO-14001 - Cert. n° 0256A/0



Tecniplast reserves the right to change or modify product and or specifications without notice or obligation.  
F-Suite is a Tecniplast trademark.