

GP-SUITE™

ENRICHED GUINEA PIG HOUSING

- **3D enriched environment:**
Optimal animal welfare
- **High quality plastics and production process:**
Easy cleaning and long lifespan
- **Available with solid or perforated bottom:**
Maximum flexibility



HOUSING



TECNIPLAST
innovation through passion™



HOUSING

GP-Suite: Enriched Guinea Pig Housing

1/Animal Welfare

- New injection moulded perforated Noryl® floor:
 - No sharp edges to **enhance animal comfort**.
 - Fast draining and drying design to maintain a **clean environment**.
 - Extra holes on the corners to **speed up the draining and drying processes**.
 - **Warm and comfortable** cage.
- A raised shelf with ramps creates a **three-dimensional enriched environment**.


2/Ergonomics and Cost-effective

- New injection moulded Noryl cage means:
 - Smooth edges that **eliminate dirty traps**, ensuring **perfect cleaning**.
 - Stacking lugs for **proper stacking**, rather than "stucking."
 - Thickness uniformity for a **lighter and more resistant cage**.
 - Stress- and tension- free plastic offering greater thermal and mechanical resistance for a **long lasting product** that withstands repetitive washing and autoclave cycles.
- Rounded corners and lighter rack ensure **easy and safe transportation**.
- Compact housing solution for an **higher density**.

3/Flexibility

Available either with solid or perforated bottom to accommodate any research need.

Animal Welfare **1**

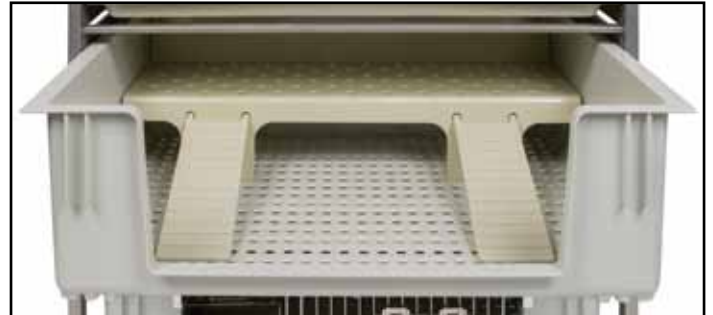


Environmental enrichment



Extra holes in the four corners

Ergonomics and Cost-effective **2**

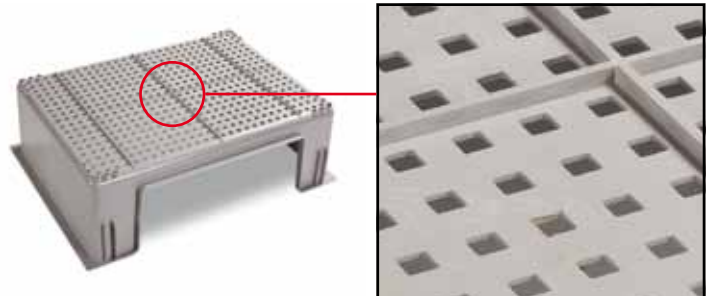

Injection moulded cage base: no dirty traps

TECHNICAL SHEET

Cage Specification and Dimensions

"W" type solid/perforated cage base

Size (L x W x H): 814 x 610 x 256 mm / 32.95 x 24.02 x 10.08 in

Floor area: 4000 cm² / 4.3 sq ft

Internal Height: 281 mm / 11.08 inch

Cage Bottom Type	Rack Configuration	Overall Rack Dimensions L x W x H - mm/in
Solid	5 cages (1 W x 5 H)	860 x 780 (846*) x 1656 33.86 x 30.71 (33.31*) x 65.20
Solid	10 cages (2 W x 5 H)	1700 x 780 (846*) x 1656 66.93 x 30.71 (33.31*) x 65.20
Perforated**	4 cages (1 W x 4 H)	860 x 780 (846*) x 1656 33.86 x 30.71 (33.31*) x 65.20
Perforated**	8 cages (2 W x 4 H)	1700 x 780 (846*) x 1656 66.93 x 30.71 (33.31*) x 65.20

* With feeders and water bottles.

** The rack with perforated cages comes complete with waste trays.



ISO-9001
Cert. n° 0875



ISO-14001
Cert. n° 0256A/0



HEADQUARTERS
ITALY | Tel. +39 0332 809711 • www.tecniplast.it • E-mail: tecnicom@tecniplast.it

SISTER COMPANIES

AUSTRALIA/NEW ZEALAND | Tel. + 61 2 8845 6500 • www.tecniplast.it • E-mail: info@tecniplast.com.au

CHINA | Tel. + 86 21 50810920 • www.tecniplast.cn • E-mail: tecniplastchina@tecniplast.it

FRANCE | Tel. 04 72 52 94 41 • www.tecniplast.fr • E-mail: info@tecniplast.fr

GERMANY | Tel. 08805 921320 • www.tecniplast.de • E-mail: info@tecniplast.de

UNITED KINGDOM | Tel. 0845 0504556 • www.tecniplast.it • E-mail: info@tecniplastuk.com

JAPAN | Tel. + 03 5770 5375 • www.tecniplast.it • E-mail: info@tecniplastjapan.co.jp

USA/CANADA | Toll Free: 877.669.2243 • www.tecniplastusa.com • E-mail: info@tecniplastusa.com

To find your local distributor, please visit www.tecniplast.it

Tecniplast reserves the right to change or modify product and or specifications without notice or obligation.
GP-Suite is a Tecniplast trademark.