

TECNIPLAST



BC-PLAST
PLASTIC BOTTLE CRATE




LOGISTICS

BC-PLAST


PLASTIC BOTTLE CRATE

Light. Robust. Thermal Resistant.

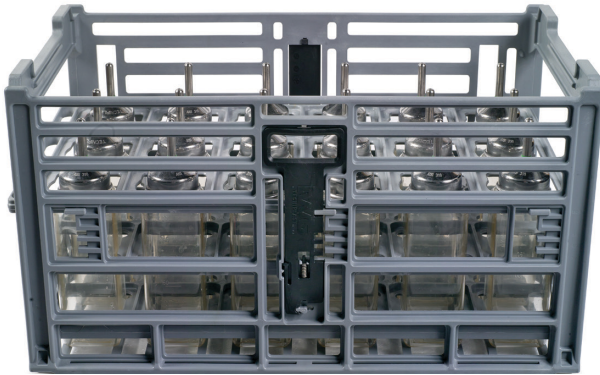

STANDARD FEATURES

- Weights less than 6.5 lbs (3 kg). Autoclavable.
- Fits bottle sizes from 13 oz to 26 oz (200 ml to 400 ml)
- Made of glass-filled injection molding PPA for structural endurance and thermal resistance.
- Embedded lid with quick lock-and-release system secures the bottles during upside down handling.
- Quick release sliding lid combined with a simple pull up mechanism provide greater ergonomics.

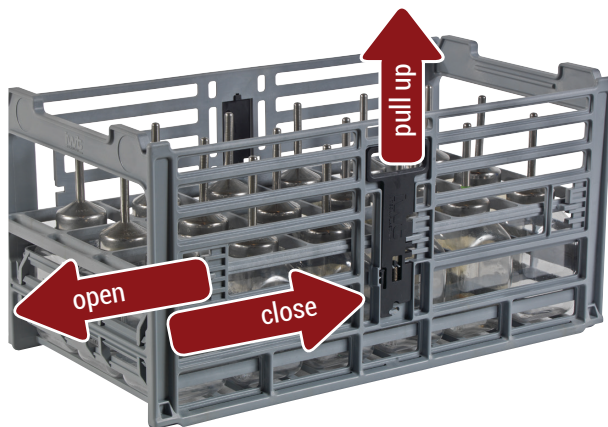

BC-Plast features an embedded lid with sliding open/close mechanism and quick lock & release system.

LIGHT IS THE NEW STRONG FOR WASHING OPERATIONS

BENEFITS

- 35% lighter than stainless steel.
- High thermal, moisture and chemical resistant.
- Compatible with most washing equipment.
- Robust structural design and construction.
- Suitable for manual and automated processes.
- Great ergonomics and improved usability.

Withstands long periods of **stress and strain** during the **day-to-day** washroom operation and autoclaving.

BC-Plast's performance and materials characteristics are proven to remain unchanged throughout the process*.



Compatible with most rack, cage and cabinet washers logistic carts and transport trolleys on the market.

(*) Consult the User's Technical Data Sheet for details about how to use this product.