


TECNIPLAST



ARIA CS48
CHANGING STATION

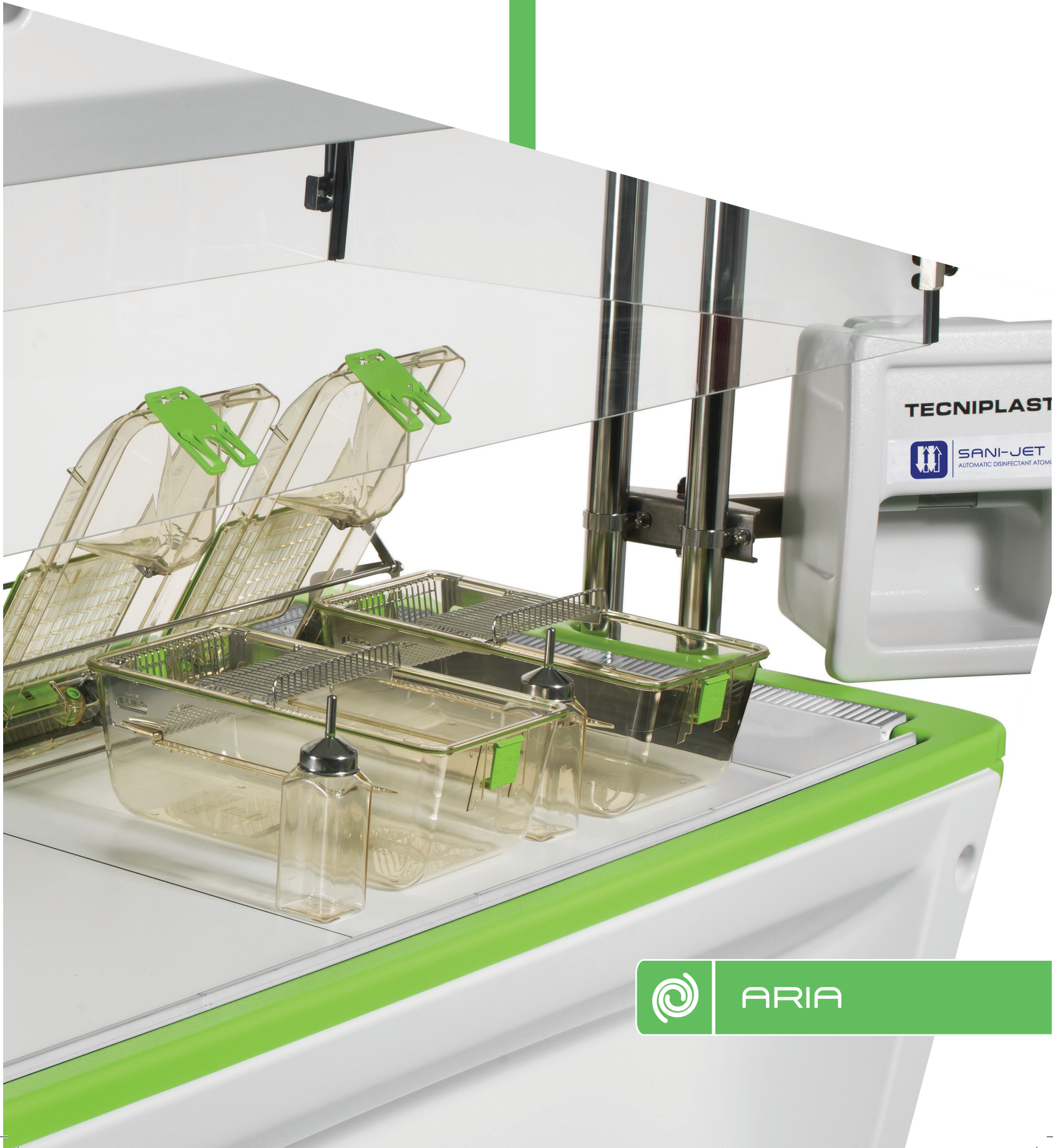

ARIA

ARIA CS48 CHANGING STATION

Dual Access. ISO Class 4 Air Quality. Innovative **Flat or Recessed** Working Surface

- Effective barrier; prevents cross-contamination and allergen exposure
- Flexibility. Dual configuration fits from the smallest to the largest cages
- Lighter. 20% less weight than other models on the market

STANDARD FEATURES

- Large working surface of size 4ft (1.2m). Fits rat and mouse cages.
- ISO Class 4 Air Quality (ISO 14644-1) for maximum protection
- Dual configuration of the working surface: flat or recessed
- Weights 407lbs (185kg); 20% lighter than others on the market
- 5" (125mm) casters for quiet and easy mobility
- Ergonomic design and sturdy construction

ENERGY SAVING FEATURES

- Energy saving blowers (ECM Technology)
- Low energy LED strip provides uniform light on the work surface
- A dimmer controls light intensity to accepted levels for albino strands
- Injection molded polysulfone airflow grille around the work surface reduces noise and turbulence generated by inside air flow

OVERALL DIMENSIONS (W x D x H):

55.70 x 31.1 x 76.9/89.68 in

1414 x 790 x 1954/2278 mm

WORKING AREA DIMENSIONS (W x D x H):

44.09 x 24 x 25.2 in (flat mode) / 27.16 in (recessed mode)

1120 x 610 x 640 mm (flat mode) / 690 mm (recessed mode)



20% lighter than other models on the market

BENEFITS

- Acoustic comfort. Low levels of noise and vibration
- Overall energy savings and reduced heat load emissions
- Safety and great ergonomics
- Light weight; easy to transport around the facility
- High cleanability
- Reduced turbulence and noise inside working surface

The ARIA CS48 is independently certified (UL/CSA) and compliant with CE requirements. Quality materials are in conformity to guarantee the strictest safety standards.

FOR INFORMATION email info@tecniplastusa.com



 **TECNIPLAST**



WWW.TECNIPLASTUSA.COM

info@tecniplastusa.com | 1.877.669.2243

info@tecniplastcanada.com | 1.855.347.8718