



Class II (Type A2):
guaranteed
protection for your
staff and
animals!



Mobile Class II Biosafety Changing Station specifically designed for Germ-Free applications

Class II Microbiological Safety Cabinet specifically designed for Animal Research Laboratories, Tecniplast ARIA ISO-TECH60 P assures an effective barrier, preserving germ free health status of the animals and preventing cross contamination as well as allergens and potential pathogens exposure.

PERFORMANCE

- **Manual Decon Lock** (option) on a side for manually decontaminating the cage before moving in and handling it or **Automatic Dunk Tank** (option) with interlocked doors which assures standardization in SOPs with germ free animals and an advanced bioprotection/ biocontainment.
- **ISO Class 3 Air Quality according to ISO 14644-1**, the aeraulic performances meets the requirement of EN 12469:2000 regulation and compliant with NSF49 as Class II Type A2
- **Highest degree of protection**: optional sliding gloves system installed along the front aperture in order to minimize the risk of human error and potential contamination
- **Transfer door (option) with or without flange** in order to introduce the autoclaved stuff or germ-free animals in complete safety.

BENEFITS

- **Very large work area** for housing all the staff needed to work with very sensitive mice
- **Enhanced ergonomics**: the improved opening height and the use of a smart sliding gloves system make any kind of manipulation simpler, more effective and safer
- **Energy-efficient**: LED light bars with dimmer and EC motor blower
- **Easy to clean**: the AISI 316 working area with removable TRESPA® work surface panels are highly resistant when carrying out routine cleaning and/or autoclaving sterilization procedures
- **Customer configurable**: different options available according to the user needs and SOPs.

STANDARD FEATURES

- **Hinged and sloped safety glass front window** allowing good visibility of the work surface and inside the cages by reducing glare
- **Electric auto lift system** (optional) for ergonomic height positioning and easy transport
- **User friendly 7" TFT Colour Touch panel** with integrated microprocessor system for easy control and operation
- **LED light bars with dimmer** (particularly beneficial for procedures with light-sensitive mice) and red LED light to work at night.

OVERALL DIMENSIONS (without exhaust box) WxDxH

with MDL and transfer door with flange:

2340 x 859 x 2007 - 2407 mm - 92" x 33"7/8 x 79" - 94"3/4

with ADT and transfer door with flange:


2215 x 859 x 2007 - 2407 mm - 87"1/4 x 33"7/8 x 79" - 94"3/4

SAFE WORKING AREA DIMENSIONS WxDxH

1495 x 490 x avg 690 mm - 58"7/8 x 19"2/8 x avg 27"1/8



MANUAL DECON LOCK

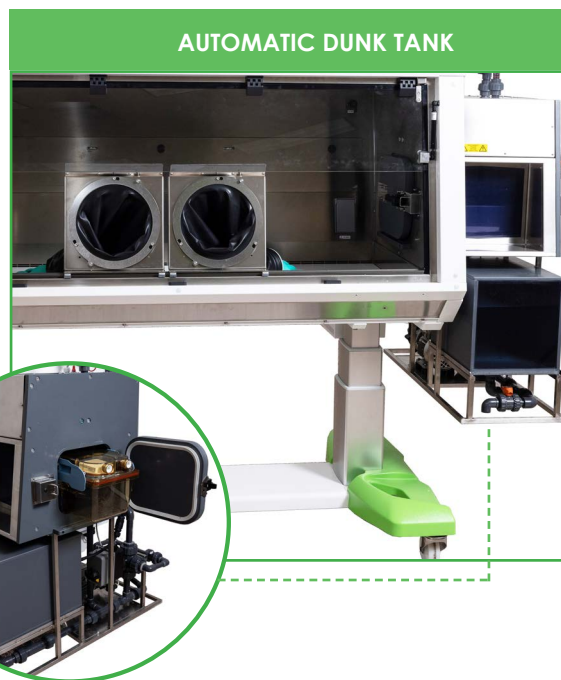

AUTOMATIC DUNK TANK



COMPANY WITH
QUALITY SYSTEM
CERTIFIED BY DNV GL
= ISO 9001 =

COMPANY WITH
ENVIRONMENTAL SYSTEM
CERTIFIED BY DNV GL
= ISO 14001 =

TECNIPLAST
innovation through passion

HEADQUARTERS

ITALY | Tel. +39 0332 809711 · www.tecniplast.it · E-mail: tecnicom@tecniplast.it; venditeitalia@tecniplast.it

SISTER COMPANIES

AUSTRALIA/NEW ZEALAND | Tel. + 61 2 8845 6500 · www.tecniplast.it · E-mail: info@tecniplast.com.au

CHINA | Tel. + 86 (0)21 50810920 · www.tecniplast.cn · E-mail: tecniplastchina@tecniplast.it

FRANCE | Tel. +33 (0)4 72 52 94 41 · www.tecniplast.fr · E-mail: info@tecniplast.fr

GERMANY | Tel. +49 (0)8805 921320 · www.tecniplast.de · E-mail: info@tecniplast.de

UNITED KINGDOM | Tel. + 44 (0)345 0504556 · www.tecniplastuk.com · E-mail: info@tecniplastuk.com

JAPAN | Tel. +81 (0)3 5770 5375 · www.tecniplast.it · E-mail: info@tecniplastjapan.co.jp

SWITZERLAND | Tel. +41 (0)81 750 31 40 · www.tecniplast.ch · Email: info@tecniplast.ch

USA | Toll Free: +1 877.669.2243 · www.tecniplastusa.com · E-mail: info@tecniplastusa.com

CANADA | Toll Free: +1 855.347.8718 · www.tecniplastcanada.com · E-mail: info@tecniplastcanada.com

CAT 208/ING - Rev Dec.2022

Tecniplast reserves the right to change or modify product and/or specifications without notice or obligation.