



D-CHAMBER


Solution for Vaporized Hydrogen Peroxide (VHP) decontamination, and Aerosolized Hydrogen Peroxide (AHP) disinfection. Also configurable for high containment (BSL3) facilities.

KEY FEATURES

- **Comfort** - Large volume chamber (6 m³ – 212 ft³) in a reduced footprint (4 m² – 43 ft²)
- **Speed** - Setup including compact Pass-Box for quick (**less than 1 hour**) decontamination of miscellaneous items
- **Easy Installation** - Minimized building requirements:
 - ▶ Pit-less design without the need for any ramp;
 - ▶ Fully serviceable either from the inside or from the dirty side for a real wall-to-wall capability (e.g., corridors).
- **Safety** - Aeration flexibility to break down H₂O₂ levels inside the chamber:
 - ▶ Double ceiling-mounted catalyst and HEPA filters without the need for any external ductwork;
 - ▶ Exhaust fan connected to the building exhaust system;
 - ▶ Building HVAC system.



- Compatible with IWT on-board dBox Plus VHP generator or with standard external generators available in the market
- Enhanced decontamination process control and reduced cycle time thanks to the dehumidifying unit integrated into the dBox Plus generator
- Built-in atomizer for fast turn-key aerosol disinfection with low-concentrated H₂O₂

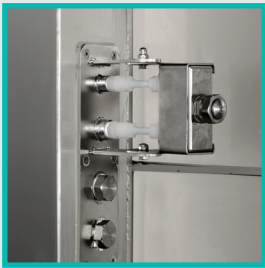

- Internal socket for power supply of the load (e.g., changing stations ventilation)
- Inclusive of low voltage distribution fan for uniform diffusion of the decontaminant during cycle
- Built-in safety:
 - ▶ Internal emergency buttons with self-releasing doors;
 - ▶ Manual sample port for ppm concentration detection prior to operator access or optional automated safety sensor for low H₂O₂ concentration inside the chamber to allow door opening under safe conditions.



COMPANY WITH
QUALITY SYSTEM
CERTIFIED BY DNV
ISO 9001

COMPANY WITH
ENVIRONMENTAL SYST
CERTIFIED BY DNV
ISO 14001

HEADQUARTERS IWT s.r.l. | Tel. +39 0332.96701 · www.iwtsrl.it · E-mail: info@iwtsrl.it

Sister Companies

AUSTRALIA/NEW ZEALAND | Tel. + 61 2 8845 6500 · www.tecniplast.it · E-mail: info@tecniplast.com.au

CANADA | Toll Free: +1 855.347.8718 · www.tecniplastcanada.com · E-mail: info@tecniplastcanada.com

CHINA | Tel. + 86 (0)21 50810920 · www.tecniplast.cn · E-mail: tecniplastchina@tecniplast.it

FRANCE | Tél. +33 (0)4 72 52 94 41 · www.tecniplast.fr · E-mail: info@tecniplast.fr

GERMANY | Tel. +49 (0)8805 921320 · www.tecniplast.de · E-mail: info@tecniplast.de

JAPAN | Tel. +81 (0)3 5770 5375 · www.tecniplast.it · E-mail: info@tecniplastjapan.co.jp

SWITZERLAND | Tel +41 (0)81 750 31 40 · www.tecniplast.ch · Email: info@tecniplast.ch

UNITED KINGDOM | Tel. + 44 (0)345 0504556 · www.tecniplastuk.com · E-mail: info@tecniplastuk.com

USA | Toll Free: +1 877.669.2243 · www.tecniplastusa.com · E-mail: info@tecniplastusa.com

CAT MKTBRC0001330 - Rev Oct2023

IWT reserves the right to change or modify product and/or specifications without notice or obligation.