




EOLUS CAGE INVERTER

EOLUS CAGE INVERTER (C.I.) IS THE PERFECT MATCH WITH YOUR IWT ARCADIA TUNNEL WASHER TO COMPLETE MODERN CAGE WASH AREAS OPERATIONAL FLOWS, FROM DIRTY CAGES TO CLEAN REFILLED ONES.



EFFICIENT PROCESSING

- **High throughput:** Patent Pending tilting bridge solution to properly turn upside down cage bases presenting them for dispensing. Up to 960 type ILL cages/hour
- **Smooth cage processing:** conveyORIZED full length flat belt to allow the transit of IVC components (Tops and Wire Lids) from TW to operator position at the end of the line
- **Dispensing accuracy:** automatic detection of cages quantity for punctual and accurate bedding dispensing

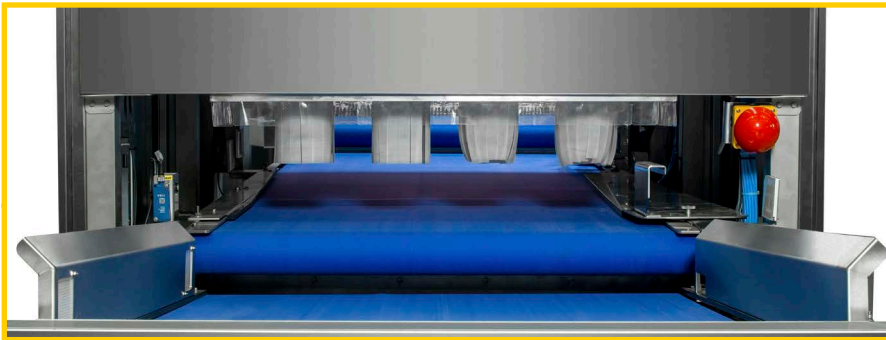


SAFETY & FLEXIBILITY

- **Operators' protection:** thanks to laser scanner technology, users can move parts without the need of physical barriers
- **Versatility:** dispensing technique via mixer/gravity technologies or mixer/augers for either free-flowing or light material processing
- **Refilling flexibility:** bedding dispenser can be refilled via centralized vacuum system or standalone compact vacuum option
- **Extra system autonomy** thanks to optional conveyORIZED buffer (1.2m – 4ft)

ERGONOMICS & EASE OF USE

- **Easy and quick items retrieval:** full frontal access to tunnel washer unloading module for direct retrieve of miscellaneous items when tilting bridge is in resting position
- **User-friendly:** operators always work at the best ergonomic height when collecting cages (800/900mm – 2 ¾/3 ft based on cage height)
- **Easy and quick cleaning procedures:** direct access to the waste tray for bedding spillage at the end of the line and fast cleaning operations via wet towel over running belt.

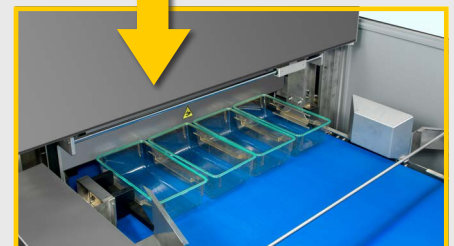

COMPANY WITH
QUALITY SYSTEM
CERTIFIED BY DNV
ISO 9001

COMPANY WITH
ENVIRONMENTAL SYST
CERTIFIED BY DNV
ISO 14001

HEADQUARTERS IWT s.r.l. | Tel. +39 0332.96701 · www.iwt srl.it · E-mail: info@iwt srl.it

Sister Companies

AUSTRALIA/NEW ZEALAND | Tel. + 61 2 8845 6500 · www.tecniplast.it · E-mail: info@tecniplast.com.au

CANADA | Toll Free: +1 855.347.8718 · www.tecniplastcanada.com · E-mail: info@tecniplastcanada.com

CHINA | Tel. + 86 (0)21 50810920 · www.tecniplast.cn · E-mail: tecniplastchina@tecniplast.it

FRANCE | Tél. +33 (0)4 72 52 94 41 · www.tecniplast.fr · E-mail: info@tecniplast.fr

GERMANY | Tel. +49 (0)8805 921320 · www.tecniplast.de · E-mail: info@tecniplast.de

JAPAN | Tel. +81 (0)3 5770 5375 · www.tecniplast.it · E-mail: info@tecniplastjapan.co.jp

SWITZERLAND | Tel +41 (0)81 750 31 40 · www.tecniplast.ch · Email: info@tecniplast.ch

UNITED KINGDOM | Tel. + 44 (0)345 0504556 · www.tecniplastuk.com · E-mail: info@tecniplastuk.com

USA | Toll Free: +1 877.669.2243 · www.tecniplastusa.com · E-mail: info@tecniplastusa.com

CAT MKTBRC000117L Rev OCT2023

IWT reserves the right to change or modify product and/or specifications without notice or obligation.