


 **EASY 400**
C A B I N E T W A S H E R

- Double loading level for extra washing capacity
- Compact size with built-in flexibility
- IWT quality at very affordable price
- Easy to install and very simple to use

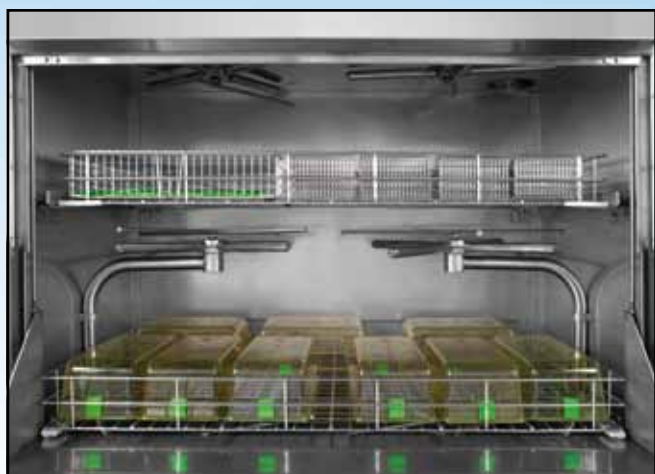

WASHING

1/ Double loading level for extra washing capacity and flexibility

EASY 400 Series Cabinet Washer is yet another product in the arsenal of engineered washing solutions by IWT.

The machine has been designed with the goal of having a good ratio footprint/load capacity.

The result is a double deck where different items, such as cages, bottle baskets, miscellaneous crates can be loaded, making this machine ideal for all that facilities where flexibility, i.e. wide loads variety, is a need.

2/IWT quality at affordable prices

EASY 400 Series Cabinet Washer features purpose-designed rotary arms on three levels to guarantee a good coverage of both decks; separate wash and rinse circuit; two wash pumps; one rinse pump and two detergent dosing pumps.

The machine is manufactured in AISI 304 stainless steel, Scotch Brite finishing.

All this comes with a very competitive price.

3/Easy to install and very simple to use

The machine comes in one piece and needs only water and electrical connections.

The possibility of placing the drain both on the floor or on the wall and the absence of an exhaust ductwork make the installation smooth and fast.

The getting started is as fast as the installation: once the items are loaded, the operator closes the door and selects the desired washing cycle just pushing a button.

Accessories

- S/s crate for 18 bottles (9LBCBM18---)
- S/s universal mesh crate for caps (9LBCBMMAGIC)
- S/s universal mesh crate with hinged lid for miscellaneous accessories (9LBCA003)
- S/s bottle crates transport trolley with drawers (ACCACBO---)

TECHNICAL SHEET

Specifications and Dimensions

Easy 400 Series Cabinet Washer					
Code	Doors	Heating System	Overall Dimensions (W x D x H) mm	Usable Wash Dimensions (W x D x H) mm	Loading Capacity (example)
9LAVR401EZ	1	Electric	1440 x 1275 x 2285	Bottom level: 1290 x 665 x 315 Upper level: 1290 x 640 x 315	18 type II cages or 6 x 18 bottle crates and 2 cap crates

Electrical supply: 400v 50Hz three-phase + N + E

Water: 3/4 °C - 2,5 to 4 bar - Max hardness CaCO₃: 150ppm (15°F) - Suggested supplied water temperature 50-60°C

Notes: Height and Depth are calculated with open sliding door.

IWT will provide mechanical information, service requirements and detailed installation drawings on request.



WASHING

Cabinet Washer Easy 400 Series

Double loading level for extra washing capacity and flexibility

1



IWT quality at affordable prices

2



Easy to install and very simple to use

3



iwt

a **TECNIPLAST** company

HEADQUARTERS **IWT s.r.l.** | Tel. +39 0332.96701 • www.iwtsrl.it • E-mail: info@iwtsrl.it

SISTER COMPANIES

AUSTRALIA/NEW ZEALAND | Tel. + 61 2 8845 6500 • www.tecniplast.it • E-mail: info@tecniplast.com.au

CHINA | Tel. + 86 21 50810920 • www.tecniplast.cn • E-mail: tecniplastchina@tecniplast.it

FRANCE | Tél. 04 72 52 94 41 • www.tecniplast.fr • E-mail: info@tecniplast.fr

GERMANY | Tel. 08805 921320 • www.tecniplast.de • E-mail: info@tecniplast.de

UNITED KINGDOM | Tel. 0845 0504556 • www.tecniplast.it • E-mail: info@tecniplastuk.com

JAPAN | Tel. + 03 5770 5375 • www.tecniplast.it • E-mail: info@tecniplastjapan.co.jp

USA/CANADA | Toll Free: 877.669.2243 • www.tecniplastusa.com • E-mail: info@tecniplastusa.com

To find your local distributor, please visit www.tecniplast.it



ISO-9001:2000 - Cert. n° 3406/1

IWT SRL Società unipersonale • Società soggetta a Direzione e Coordinamento di Tecniplast S.p.A.
Sede legale e amministrativa: via Galliani 68b 21020 CASALE LITTA (VA) ITALY • C.F. e P.IVA: IT 02002720122

W-DP052/ING - Rev. February 2014