



# NHP HOUSING AND WASHING

## Marmoset Housing

European Standard

Built with welfare and flexibility in mind.  
Fully washable.

### VERTICAL HEIGHT AND BUILT PROMOTE NATURAL ARBOREAL BEHAVIOR

The Tecniplast Marmoset Family cage offers the greatest vertical height PLUS other features that facilitate and promote the animals' natural arboreal behavior.

- Height on the inside is 79" (6.58 ft.)
- Stainless steel mesh in all six sides
- True nest box located at the top of the cage

Fig. 1




79"

### TRUE NEST BOX - MULTIPURPOSE

A nest box made of plastic is located at the top of the cage (Fig. 1). This box:

- Provides protected shelter for the animals
- Promotes their natural climbing behavior
- Facilitates their capture and transport

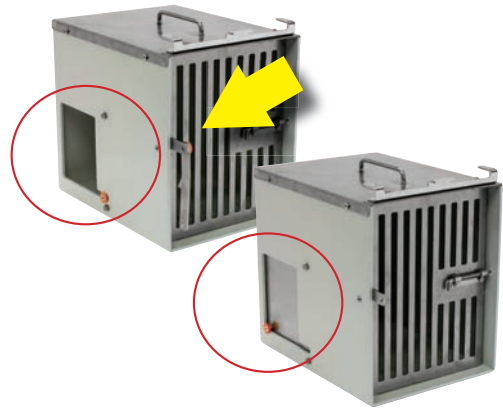

### EASY CAPTURE AND TRANSPORT

can be done by one person. The latch system and design allow the operator to:

- Open and close box's gate from the outside without opening the cage's front door (Fig. 1)
- Easily slide the box out of the cage without using any tools.
- Safely transport the animals.

## FULLY WASHABLE

- Marmoset cages fit right in the ATLANTIS XL rack washer chamber.
- The ATLANTIS XL is a high-efficiency rack washer with customizable washing cycles.
- Trays and accessories can be washed using a dedicated presentation rack that goes directly into the chamber.

## SCALABILITY MEETS THE NEEDS FOR INDIVIDUAL AND GROUP HOUSING

- Securely connect an unlimited number of cages to house large groups.
- Convert from individual housing for study conditions to family or breeding sized units by utilizing dedicated plastic sliding panels.

## ENRICHMENT - THE CAGE DESIGN FACILITATES ALL COMMON ENRICHMENT

- Mesh in all six sides promotes the animals' natural climbing behavior.
- Vertical height and built-in facilitate enrichment.
- Nest box located at the top promotes the mothers' natural arboreal behavior.



# ATLANTIS HIGH EFFICIENCY CAGE WASHER



## SLIDING PANELS ADD FLEXIBILITY FOR GROUP BEHAVIOR STUDIES

- Block visual contact between groups.
- Divide a cage into two full-height cages preserving the cage's climbing attribute and inside dimensions.
- Install or remove without disturbing the animals.
- Toll free removal or install by only one operator.

## OTHER OPTIONS AVAILABLE

### CUSTOMIZED

Group playgrounds available upon request.



## ATLANTIS XL FOR NHP CAGES, PARTS, AND ACCESSORIES

### FEATURES:

- High-efficiency NHP cage washer
- 20-minute washing cycle.
- Reduced utility consumptions.
- Modular design and less infrastructure requirements for easy installation.
- Hot air drying with heat-recovery system. Automatic tilting floor.
- Double washing tank for alkaline and acid phases.
- H2O2 decontamination. AWS flushing and steam thermal disinfection.
- Optional configurations for energy recovery.
- Oscillating arms and vertical manifold to achieve full coverage.

### COMPARED TO OTHER WASHERS IN THE MARKET, THE ATLANTIS:

1. Delivers 50% more water flow and pressure for effective cleaning.
2. Requires 10% to 20% less footprint, including maintenance access space.

ATLANTIS XL is fully certified as compliant with AAALAC requirements.