



NHP HOUSING AND WASHING

Macaque Housing

Cage Group 3 & Cage Group 4

Built with welfare and flexibility in mind.
Fully washable.

MODULAR INTERLINKING to expand

FROM INDIVIDUAL TO GROUP HOUSING

- Interlinking capability with secure locking system.
- Removable sliding panels made of Trespa® add comfort for the animals.
- Panels lighter than steel. Easy to manipulate.
- Great durability. Easy to clean.

STANDARD FEATURES

- Dual cage configuration
- Interlinking capability
- Waste tray below each cage
- Squeeze back for restraint
- Protected bottle holder
- Easy to clean feeder
- Horizontal sliding mesh door
- Removable Trespa panels
- S/s back panel




Secured multi-linking capability



Cage Group 3

Floor Area: 4.3 ft² (0.4 m²)

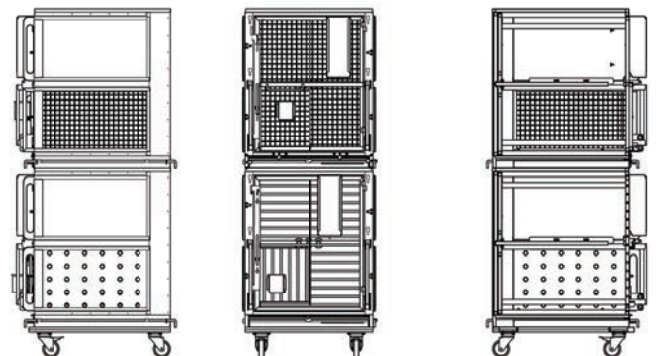

FULLY WASHABLE

- Macaque cages fit right into the ATLANTIS XL rack washer chamber.
- The ATLANTIS XL is a high-efficiency rack washer with customizable washing cycles.
- Trays and accessories can be washed using a dedicated presentation rack that goes directly into the chamber.

INNOVATIVE BOTTLE HOLDER

Encased water bottle holder that keeps the bottle protected from being removed or destroyed by the animals.

FLEXIBLE BUILT

ADD OR MODIFY FEATURES SUCH AS:

- Auto-watering
- Forest or plain Trespa panels
- Sliding or hinged doors
- Solid or perforated panels
- Transfer box

ACCESSORIES

- Animal transfer cart
- Pan cart for transport
- Pan presentation rack
- Enrichment




Encased water bottle holder

CUSTOMIZED

Group playgrounds available upon request.

Cage Group 4

Floor Area: 6.0 ft² (0.56 m²)

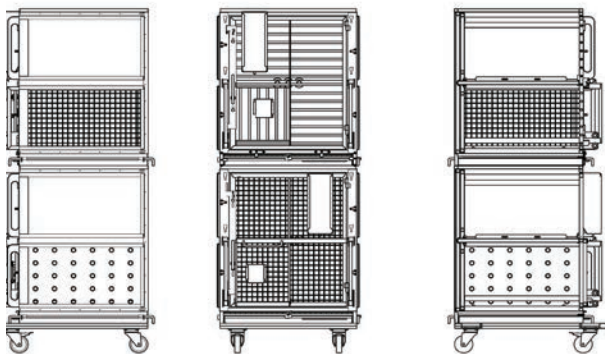

ATLANTIS
HIGH
EFFICIENCY
CAGE
WASHER

NHP HOUSING AND WASHING

Marmoset Housing

European Standard

Built with welfare and flexibility in mind.
Fully washable.

VERTICAL HEIGHT AND BUILT PROMOTE NATURAL ARBOREAL BEHAVIOR

The Tecniplast Marmoset Family cage offers the greatest vertical height PLUS other features that facilitate and promote the animals' natural arboreal behavior.

- Height on the inside is 79" (6.58 ft.)
- Stainless steel mesh in all six sides
- True nest box located at the top of the cage

Fig. 1




79"

TRUE NEST BOX - MULTIPURPOSE

A nest box made of plastic is located at the top of the cage (Fig. 1). This box:

- Provides protected shelter for the animals
- Promotes their natural climbing behavior
- Facilitates their capture and transport

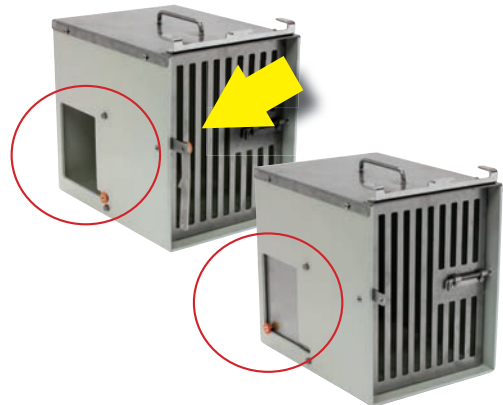

EASY CAPTURE AND TRANSPORT

can be done by one person. The latch system and design allow the operator to:

- Open and close box's gate from the outside without opening the cage's front door (Fig. 1)
- Easily slide the box out of the cage without using any tools.
- Safely transport the animals.

FULLY WASHABLE

- Marmoset cages fit right in the ATLANTIS XL rack washer chamber.
- The ATLANTIS XL is a high-efficiency rack washer with customizable washing cycles.
- Trays and accessories can be washed using a dedicated presentation rack that goes directly into the chamber.

SCALABILITY MEETS THE NEEDS FOR INDIVIDUAL AND GROUP HOUSING

- Securely connect an unlimited number of cages to house large groups.
- Convert from individual housing for study conditions to family or breeding sized units by utilizing dedicated plastic sliding panels.

ENRICHMENT - THE CAGE DESIGN FACILITATES ALL COMMON ENRICHMENT

- Mesh in all six sides promotes the animals' natural climbing behavior.
- Vertical height and built-in facilitate enrichment.
- Nest box located at the top promotes the mothers' natural arboreal behavior.



ATLANTIS HIGH EFFICIENCY CAGE WASHER



SLIDING PANELS ADD FLEXIBILITY FOR GROUP BEHAVIOR STUDIES

- Block visual contact between groups.
- Divide a cage into two full-height cages preserving the cage's climbing attribute and inside dimensions.
- Install or remove without disturbing the animals.
- Toll free removal or install by only one operator.

OTHER OPTIONS AVAILABLE

CUSTOMIZED

Group playgrounds available upon request.



ATLANTIS XL FOR NHP CAGES, PARTS, AND ACCESSORIES

FEATURES:

- High-efficiency NHP cage washer
- 20-minute washing cycle.
- Reduced utility consumptions.
- Modular design and less infrastructure requirements for easy installation.
- Hot air drying with heat-recovery system. Automatic tilting floor.
- Double washing tank for alkaline and acid phases.
- H2O2 decontamination. AWS flushing and steam thermal disinfection.
- Optional configurations for energy recovery.
- Oscillating arms and vertical manifold to achieve full coverage.

COMPARED TO OTHER WASHERS IN THE MARKET, THE ATLANTIS:

1. Delivers 50% more water flow and pressure for effective cleaning.
2. Requires 10% to 20% less footprint, including maintenance access space.

ATLANTIS XL is fully certified as compliant with AAALAC requirements.