

TECNIPLAST



NHP HOUSING
MACAQUES



HOUSING

NHP - Macaques Group 3 and 4 Cages

The Tecniplast NHP Macaque Group 3 and 4 housing meet the requirements in the NIH Guide, and is built with welfare and flexibility in mind.

MODULAR INTERLINKING to expand FROM INDIVIDUAL TO GROUP HOUSING

- Interlinking capability with secure locking system
- Removable sliding panels made of Trespa® add comfort for the animals.
- Panels lighter than steel. Easy to manipulate.
- Great durability. Easy to clean.



INNOVATIVE BOTTLE HOLDER

As a unique feature, the Tecniplast cage features an encased water bottle holder that keeps the bottle protected from being removed or destroyed by the animals.

FLEXIBILITY to meet your needs


To your convenience, add or modify features such as: Auto-watering, forest or plain Trespa panels, solid or perforated panels, and a transfer box.

Standard Features

- Dual cage configuration
- Interlinking capability
- Waste tray below each cage (fig.1)
- Squeeze back for restraint
- Protected bottle holder
- Easy to clean feeder
- Horizontal sliding mesh door
- Removable Trespa panels
- S/s back panel



Fig.1

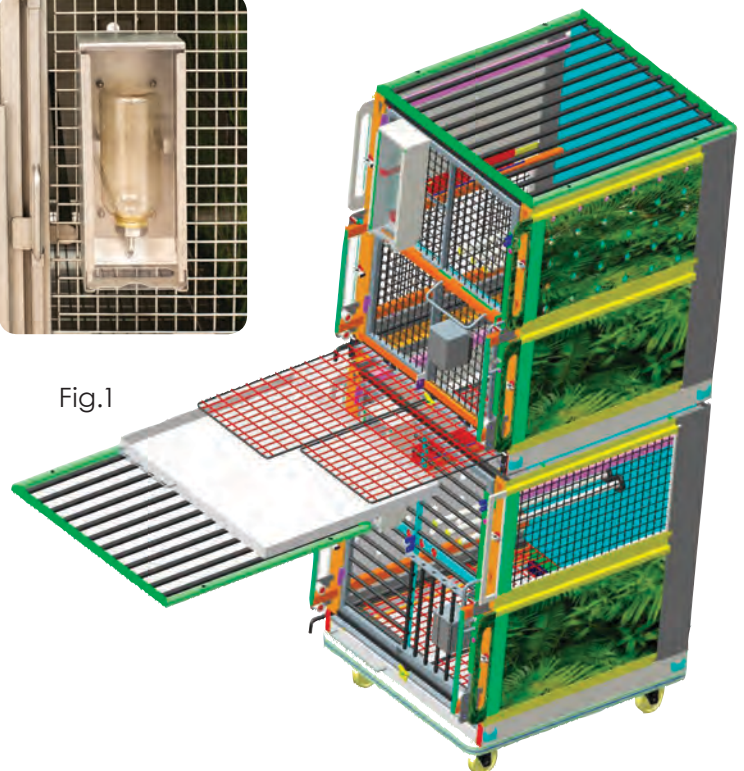



Options


- One or two front doors, sliding or hinged; bars or mesh
- Top panel: Mesh or bars
- Back Panel: Trespa
- Side Panels:
 - Solid white Trespa or forest look alike
 - Perforated White Trespa
 - Perforated Forest Trespa
 - S/s Mesh
- Auto-watering
- One tray per cage
- Customized group playground area (fig.2)

Accessories

- Animal transfer cart
- Pan cart for transport
- Pan presentation rack
- Enrichment




Fig.2

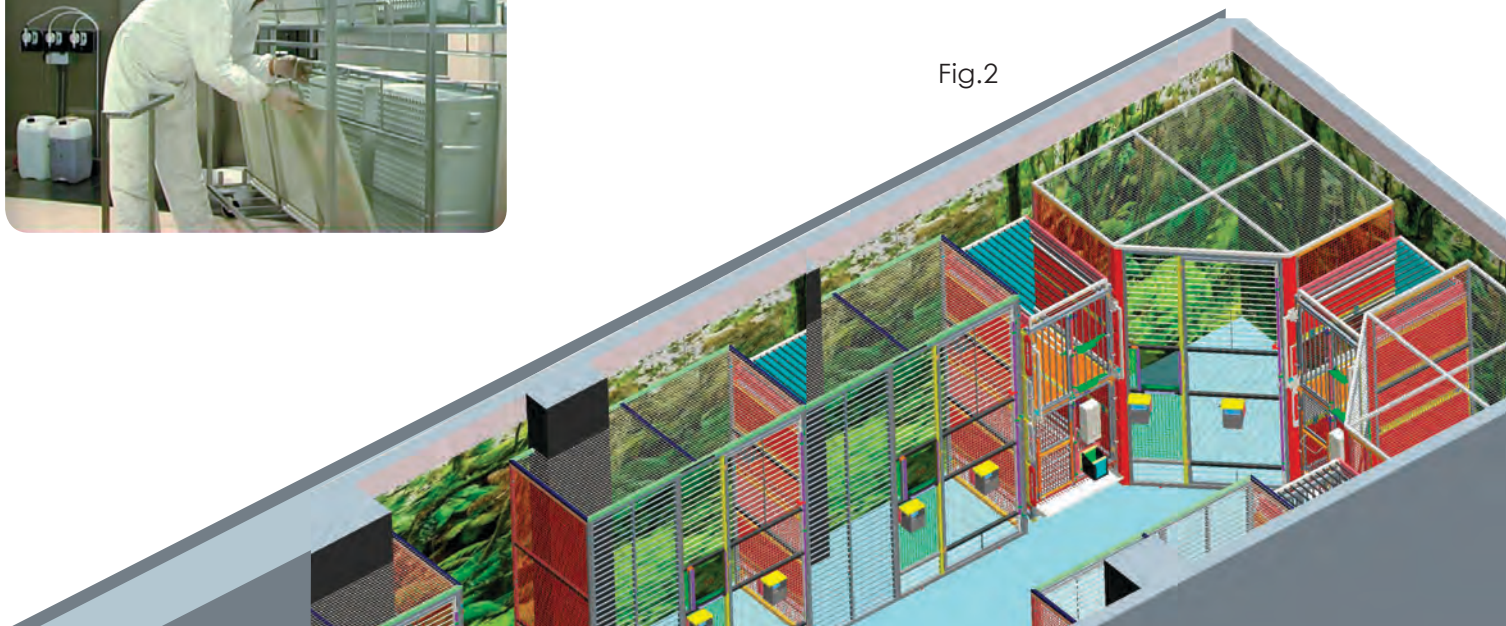


Cage Group 3

Floor Area: 4.3 ft² (0.4 m²)

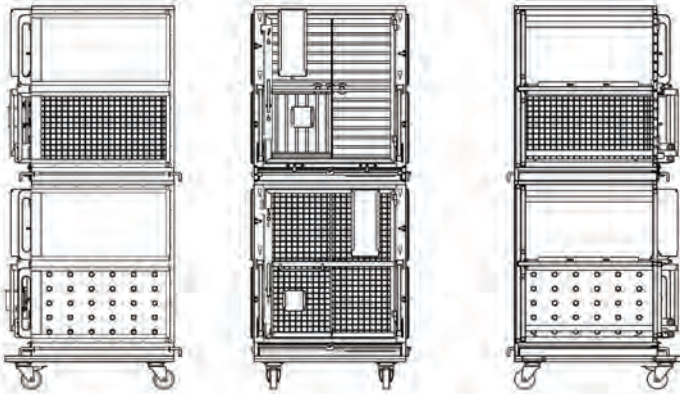

Cage Group 4

Floor Area: 6.0 ft² (0.56 m²)

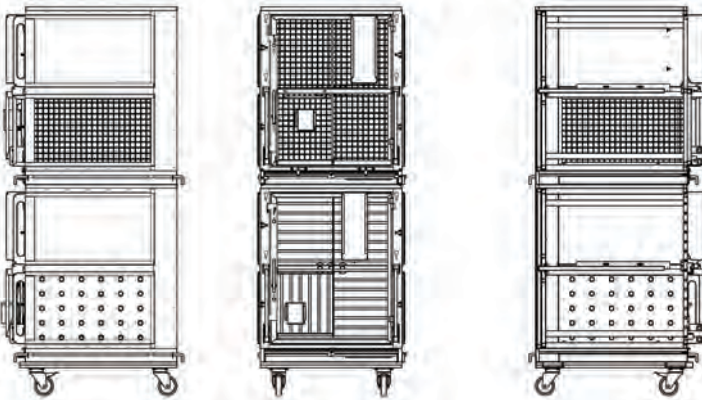

Also available NHP Marmoset Family Housing

- Height on the inside is 79 inches (6.58 ft.)
- Stainless steel mesh in all six sides
- Facilitates all common enrichment
- True nest box located at the top of the cage
- Removable side panels
- Plastic sliding dividers

Ask for details



NHP_Macaque_July 2019

 **TECNIPLAST**

COMPANY WITH
QUALITY SYSTEM
CERTIFIED BY DNV GL
= ISO 9001 =

COMPANY WITH
QUALITY SYSTEM
CERTIFIED BY DNV GL
= ISO 14001 =

USA: www.tecniplastusa.com
info@tecniplastusa.com | 1.877.669.2243

CANADA: www.tecniplastcanada.com
info@tecniplastcanada.com | 1.855.347.8718